


Vilim Ball (1.1)

vibration physiotherapy device



VILIM ball

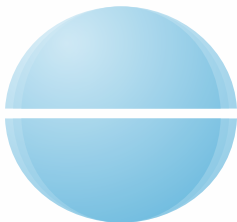









Instructions for Use

Rev6

Contents

1. Meaning of symbols	3
2. Introduction	4
3. Intended application of the device	4
4. Contraindications and precautions	4
5. Serial number description	6
6. Device description	6
7. Preparation for the use	7
8. Operating the device	8
9. Personalized therapy	11
10. Safety requirements	12
11. Package contents	13
12. Recommendations for application	14
13. Transportation and storage	14
14. Disposal of the device	15
15. Cleaning	16
16. Troubleshooting	17
17. Guarantee of the manufacturer	18
18. Technical specifications	18

1. MEANING OF SYMBOLS

	WARNING: You must read these warnings before using Vilim Ball.
	CAUTION: Contains important information about the operation and maintenance of the Vilim Ball. Read it carefully in order to avoid any problems.
	RECOMMENDATION: Contains recommended information for help with operation of the Vilim Ball.
	NOTICE: Take a user notice.
	Manufacturer
	Serial number
	Type BF applied part
	WEEE – (Waste Electrical and Electronic Equipment)
	Consult instructions for use
	CE marking of conformity (Product conforms to the essential requirements of European MDD 93/42/EEC)

2. INTRODUCTION

Vilim Ball device is non-invasive portable physiotherapeutic remedy.

The device is intended to use in home or professional environment as a mechanical vibration therapy device to temporary reduce hand tremor caused by Essential tremor or Parkinsonian tremor.

Device is patented, certificated, corresponds with requirements of Europe Union.

The effectiveness and safety have been proven by scientific research and clinical evaluation.

3. INTENDED APPLICATION OF THE DEVICE

Vilim Ball non-invasive device intended for temporary alleviation and or treatment of essential and parkinsonian tremor of upper limbs.

4. CONTRAINDICATIONS AND PRECAUTIONS

Contraindications:

- Pregnancy
- Acute thrombotic process (myocardial infarction, (acute vascular constriction)
- Implants in activated regions of the body (e.g. artificial joints)
- Acute inflammation of the locomotor system active arthrosis or arthropathy e.g. acute inflammation or swelling of joints
- Acute tendinopathy in activated regions of the body (acute tendon inflammation)
- Acute desmopathy (acute problems at the intervertebral disc)
- Fresh fractures in activated regions of the body
- Post-surgery wounds and fresh wounds in activated regions of the body or incomplete wound healing
- Rheumatoid arthritis
- Epilepsy.

Precautions:

- The Vilim Ball is not intended for other uses: cranial application to treat headache (should not be placed anywhere to induce vibrations in the head area), head or neck tremor or other conditions. It is to only be used as a handheld device to produce local vibration in the upper limb.
- The device is not created to be thrown or used in other unintended applications (e.g. heating in the microwave, chewing or gnawing). It is not intended to be a children's toy.
- The device may induce or exacerbate Raynaud's syndrome. In such event device use should be discontinued. The device induces Raynaud's syndrome is a reversible condition. Consultation of a physician is to be sought in such case.
- The device may induce or exacerbate Carpal tunnel syndrome. Carpal tunnel syndrome is a condition that causes pain and numbness in the fingers and hands, and sometimes the arms. It happens when a nerve in the wrist called the "median nerve" gets pinched or squeezed. In such case the device use should be discontinued. No irreversible damage should be expected.
- Itching redness of skin. If such condition develops and seems to be related with the device use, the use should be discontinued, and medical attention sought.

5. SERIAL NUMBER DESCRIPTION

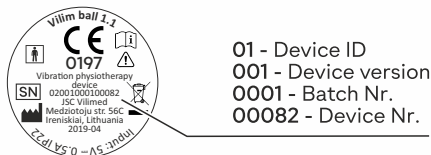


Figure 1. Device label

6. DEVICE DESCRIPTION

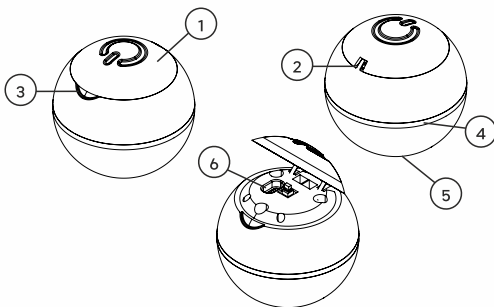


Figure 2. External view of the device

- | | |
|---------------------------|---------------------------------|
| 1. Button cap | 4. User Interface Window |
| 2. Hook for holding strap | 5. Label (bottom of the device) |
| 3. Opening cutout | 6. Charging connector |

7. PREPARATION FOR THE USE

Vilim Ball device should be held in palm and softly compressed with fingers (figure 3). Relax your hand for the best therapy results. The device can be portable / body-worn if it is turned off.

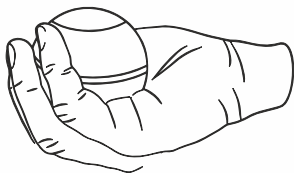








Figure 3. Holding instruction

Vibration therapy starts when Button cap (1) is pressed. Device can be stopped any time by pressing the same button cap (1). Therapy stops automatically after 10 minutes.

	Before starting the procedure, it is necessary to consult Instruction for Use of the device.
	CAUTION: Vilim Ball should not be used on patients with implanted heart stimulator / No access for persons with pacemakers.
	WARNING: Before plugging the device to charge it is necessary to check the power supply cable for mechanical damage.
	WARNING: Users with visual impairments shall have attention of caretaker for full duration of therapy.
	RECOMMENDATION: If the device was in premises with temperature lower than 5°C, then it is necessary to keep it in the workplace for at least 15 minutes before connecting to the power supply.
	RECOMMENDATION: Sounder for visually impaired can be installed to the device. For more information contact manufacturer.

8. OPERATING THE DEVICE

- The device is turned on by pressing the Button cap (1). GREEN light on User Interface Window (4) indicates that therapy is in progress.
- After 10 minutes of therapy device will stop automatically. Blue light will appear on User Interface Window (4) to inform about end of therapy.
- After the therapy Button cap (1) must be pressed to turn the device off. If Button cap (1) is not pressed device turns off automatically 5 minutes after end of therapy.
- No light in User Interface Window (4) indicated that device is turned off.

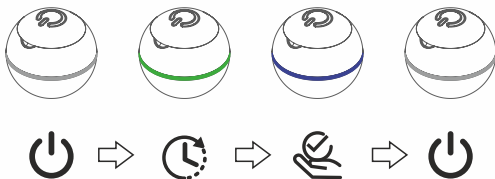


Figure 4. User Interface Window (4) instruction (off - therapy - completed - off)

- Low battery level is indicated by blinking BLUE light in User Interface Window (4) before start of therapy.
- Cap (1) must be opened in order to access the Charging connector (6). This shall be done by putting finger in the Opening cutout (3) and softly lifting Cap (1) until it opens (figure 5).

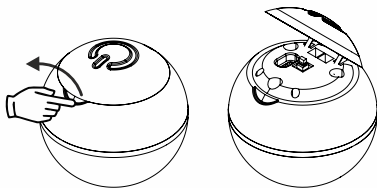


Figure 5. Opening the cap to access Charging connector (b)

- Device charging process is performed in such order:
 1. UK plug adapter (depends on region) is connected to the power supply (see figure 6).
 2. USB connector of magnetic charger cable is connected to power supply's USB socket (see figure 6).

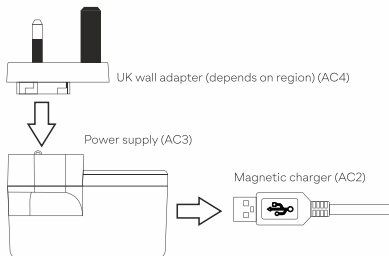


Figure 6. Preparation for charging (b)

3. Power supply is connected to standard AC outlet.
4. Magnetic connector of magnetic charger cable is connected to Vilim Ball's Charging connector (b).

- While battery is charging, BLUE light is fading in User Interface Window (4).
- If charging connector is connected and BLUE light in User Interface Window (4) turned off - charging is completed.
- If charging connector is connected and BLUE light in User Interface Window (4) is not appearing - battery is empty. In such case charging connector shall be connected for at least hour.
- Holding strap shall be fixed to a device using a simple knot on Hook for holding strap (2) as described in figure 7.

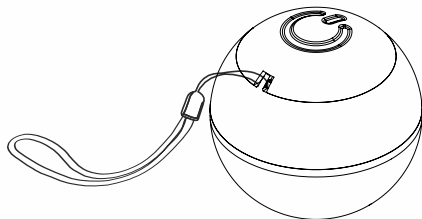




Figure 7. Holding strap.

	<p>NOTICE: Constantly blinking RED light in User Interface Window (4) means that possible error occurred. In this case, the therapy is discontinued, and the device should be switched off.</p>
	<p>WARNING: After fixing the holding trap, check if fixation is strong before use. This can be done by softly pulling the strap with one hand while holding the turned off device in other hand.</p>

9. PERSONALIZED THERAPY

Device implements system that measures individual patient's tremor. Depending on those measurements, device settings are adjusted to fit needs of each patient individually.














Take into consideration when using the Vilim Ball:





- Minimum amount of two weeks of everyday use is needed for system to fully personalize the therapy.
- Sharing device with other patients corrupts personalization data. Time to fully personalize the therapy in such case may be prolonged drastically.



NOTICE: Therapy personalization may cause sudden change of device vibration.

10. SAFETY REQUIREMENTS

	WARNING: Before plugging the device to charge it is necessary to check the power supply cable for mechanical damage.
	WARNING: It is forbidden to roll up the hook near the plug in and charger.
	WARNING: Before starting a procedure Vilim Ball device must be held on the patient's palm and only then the Button Cap can be pressed.
	WARNING: It is forbidden to use the device while battery is charging.
	WARNING: The device must be kept out of the reach of children.
	WARNING: The device is not recommended to use for children.
	WARNING: The device shall not be held above the legs to prevent damage in case of accident fall.
	WARNING: Holding strap accessory shall be used only with attention of caretaker for full therapy duration.
	WARNING: Any modification of the device is forbidden.
	WARNING: Unauthorized disassemble and any modification of device is forbidden. In case of malfunction and failure of the device always contact equipment manufacturer.
	CAUTION: Keep the device in dry area, away from sun rays.
	CAUTION: It is recommended to disinfect surface with 70 % spirit before every procedure.
	CAUTION: The accessory cables may cause cord strangulation for children.



	RECOMMENDATION: After continuous 1 hour work it is recommended to turn off the device for 15 minutes.
	RECOMMENDATION: Use device only in 5°C...+40°C ambient temperature range.
	NOTICE: Device must be periodically charged even if not used (every 3 weeks).
	NOTICE: Battery lifetime is 2 years (300 cycles), based on datasheet of the battery. The battery of the device defines lifetime of the device due to lowest lifetime of the main component. Calculations assume that an average user recharges Vilim Ball after 2-3 days of usage. That is around 147 cycles per year. Battery can be replaced after two years of usage to ensure full battery capacity.

11. PACKAGE CONTENTS

The complete device set includes:



- Vilim Ball device
- Instruction for Use
- Holding strap (Accessory marked “AC1”)
- Charger :
 - Magnetic charger cable (Accessory marked “AC2”)
 - Power supply (Accessory marked “AC3”)
 - UK wall adapter (depends on region) (Accessory marked “AC4”)
 - Power supply (AC3) specifications:

Name: UES06WNCPU-050100SPA
 Manufacturer: Dongguan Shilong Fuhua Electronics Co., Ltd
 Type: AC – DC switching power adaptor
 Input voltage: 100 Vac – 240 Vac
 Input current: 0.2A
 Input frequency: 50 Hz – 60 Hz
 Output voltage: 4.75 Vdc – 5.25 Vdc
 Output current: 1A

	CAUTION: It is forbidden to use other power supplies to charge Vilim Ball.
	CAUTION: It is forbidden to use or interconnect device with accessories, parts and materials that are not described in instructions for use.




12. RECOMMENDATIONS FOR APPLICATION

Proposed therapy time is 10 minutes per hand, three times daily, at intervals of 4 hours.

	CAUTION: Do not use device for more than 30 minutes per a day for one hand.
	WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30cm (12 inches) to any part of the medical device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

13. TRANSPORTATION AND STORAGE


Storage and transportation environments for device and its accessories must be clean and dust-free.

	CAUTION: The device can be transported in any covered transport mean, an ambient temperature range of -25 °C to + 72 °C, a relative humidity range of 10 % to 90 %.
	CAUTION: Transporting the device the stationary position must be guaranteed to avoid mechanical damage, vibration and sudden motions.
	CAUTION: For storage the device must be packed in plastic bag, placed in cardboard box and stored indoor at an ambient temperature range of + 5 °C to + 40 °C, a relative humidity range of 15 % to 90 %.

14. DISPOSAL OF THE DEVICE

Disposal of Vilim Ball non-invasive device must be performed in accordance with Directive 2012/19/EU of the European Parliament and of the Council of 4 July 2012 on waste electrical and electronic equipment (WEEE). The device must be consigned to the companies engaged in recycle and utilization of electronic waste.

- Material specifications:
 - a) Steel
 - b) Stainless Steel
 - c) Copper
 - d) Aluminum
 - e) Carbon
 - f) Bronze
 - g) Brass
 - h) Composition of polycarbonate and ABS (PC+ABS)
 - i) Silicone rubber
 - j) Polyamide (PA)
 - k) Polyurethane
- Hazardous components to be separated at the end of life cycle:
 - a) Li- ion battery
 - b) Printed Circuit Board (PCB)
 - c) Electric cables

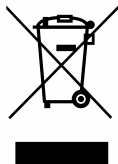



15. CLEANING

Device cleaning procedure shall be performed in this order:

- 1) Moisten a soft napkin with 1 % chloramine solution.
- 2) Take the device and close the Button cap.
- 3) Clean the device by carefully rubbing moistened napkin to its surface.
- 4) Open the Button cap.
- 5) Clean the device by carefully rubbing moistened napkin to its surface under the Button cap.
- 6) If holding strap is attached, clean it by carefully rubbing moistened napkin to its surface.
- 7) If device was being charged, take the charger and clean it by carefully rubbing moistened napkin to its surface.
- 8) Leave the device and accessories in dry place for 10 minutes.
- 9) Close the button cap.

If device and its accessories are used only by one individual, cleaning procedure shall be performed at least once in a month.

If device and its accessories are used by multiple individuals, cleaning procedures shall be performed after each use.

	WARNING: It is forbidden to disassembly the device during cleaning.
	WARNING: Device shall be kept away from dusts and cleaned immediately if contaminated.

16. TROUBLESHOOTING

Failure	Causes and actions
Battery does not charge	<ol style="list-style-type: none">1. No contact. Clean or reattach the cable.2. Battery is out of order. Contact the manufacturer.
If the green LED is not illuminating when the device is turned on (the device is loaded)	Failure of the device. Contact the manufacturer.
RED LED is constantly blinking and/or audio signal is generated	Possible internal malfunction. Contact the manufacturer.
Device works unusually noisy or vibration is not “clean”	Possible device body failure or contamination (dust, hair or other material). Contact the manufacturer.

17. GUARANTEE OF THE MANUFACTURER

A warranty period of the device and accessories in operation - 24 months from the date of realization, and a warranty period of storage - 48 months from the date of manufacturing.

The manufacturer eliminates free-of-charge all malfunctions of the device during a warranty period - 24 months at absence of mechanical damages and safeties of a seal of the manufacturer.

For assistance, if needed for using or maintaining the device or in case of the device going out of order and with suggestions or requests please contact:




JSC „Vilimed“

Vilimed.com

Medziotoju str. 56C, Ireniskiai, Lithuania (LT-53275)

Tel.: +370 646 22334

E-mail.: info@vilimed.com

	NOTICE: Service life of device and accessories is 24 months from the date of realization.
	NOTICE: Device battery needs to be replaced after 24 months from the date of realization.
	WARNING: It is strictly forbidden to repair the device yourself. In this case, the company is not liable for the consequences.

18. TECHNICAL SPECIFICATIONS

Charging input voltage:	5V DC
Charging input current:	200 mA
Battery type	Li- ion
Battery capacity:	430 mAh
Battery Voltage	3.7V
Max. battery voltage:	4.2 V
Battery compliance:	EN IEC 62133, UN 38.3
Active battery lifetime	1 – 5 hours

Stand by battery lifetime	28 days
Charging time	2 hours
Max power:	2 W
Max current:	0.7 A
Weight:	120 g
Diameter:	70mm
IP certified:	IP22*
Communication frequency:	2.4 GHz (BLE)
Medical device class:	IIa
Vibrations type:	Mechanical
Vibrations frequency range:	8Hz...18Hz

Environmental conditions for normal device operation:

Ambient temperature	5°C...+40°C
Relative humidity	15%...90%
Atmospheric pressure	700hPa...1060hPa

Environmental conditions for device transportation:

Ambient temperature	-25°C...+70°C
Relative humidity	10%...90%
Atmospheric pressure	500...1060hPa

Standards conformity:	EN IEC 60601-1
	EN IEC 62133
	EN IEC 62304
	EN 62366-1
	EN 300 328 v2.1.1
	EN ISO 10993-1
	WEEE

*Ingress protection (IP) rating describes the level of protection from solid objects and the level of protection from liquids. The Vilim Ball device is certified IP22 level protection:

Protected against solid objects over 12 mm (fingers, sticks, etc.).

Protection against direct sprays of water up to 15 degrees from the vertical position (mist, vapor, etc.).

